

CRYOPORT EXPRESS® C3™ SHIPPING SYSTEMS



Shipment Loading Instructions


NEED HELP? CONTACT CRYOPORT SYSTEMS' CLIENT CARE AT (949) 470-2305 OR EMAIL CS@CRYOPORT.COM

Step 1

Be sure to wear the appropriate personal protective equipment (PPE) when handling temperature sensitive shipments.




Step 2

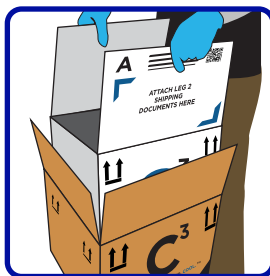

Carefully cut the top of the outer brown box and remove the inner white box by pulling the two box flaps. Remove the packing slip from the outer brown box if necessary. Throw away the outer brown box.

NOTE: If red serialized tamper-evident seals are included, record the serial numbers on all relevant documentation prior to removal. General Purpose (GP) shipping systems will arrive in a white and black box, while the Advanced Therapy Shipper (ATS™) will arrive in a white and blue box.

Step 3

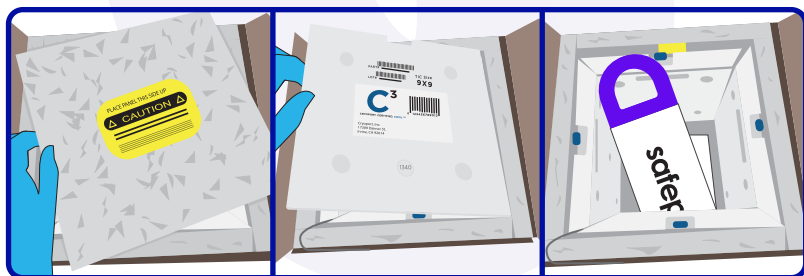

Remove the cooler lid and Phase Change Panel (PCP) to access the specimen chamber and load your commodity.

NOTE: Check the shipping system for punctures. If you see any damage, please call Cryoport Systems' Client Care at (949) 470-2305.

Step 4




Return the Phase Change Panel (PCP) to where it was before and then reinsert the cooler lid.

Step 5

For International Shipments ONLY




Remove the commercial invoices located in the shipping pouch. Then sign and date.

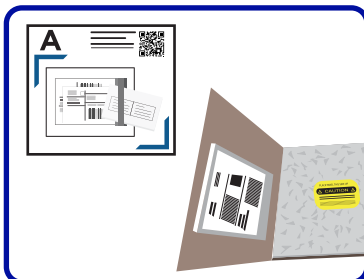

Return all shipping documents behind the Air Waybill, including any permits, forms, licenses, etc.

Step 6

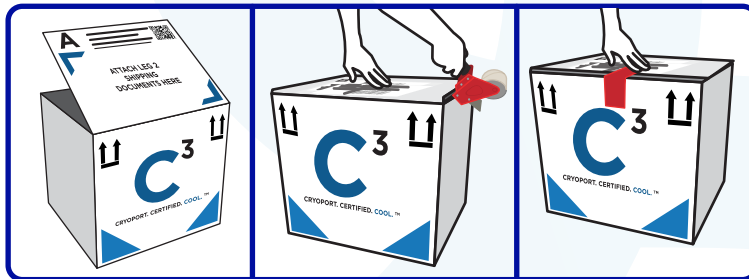

Close and seal the top of the white box with packaging tape so that flap "A" is visible on the outside of the box.

NOTE: If using the optional red serialized tamper-evident seal (ACC-9041), place it across the top of flap "A". If desired, additional seals may be ordered and applied over other openings along flap "A" for added security.