


Standard operation of the SmartPak II™ is to record every 5 minutes and to transmit all of the data at hourly intervals. The SmartPak II™ device must be in communication with a local cell tower in order to transmit data every hour.

To retrieve the most updated sensor data information, push the "Touch Sensitive Button" on the SmartPak II™ device. The most updated sensor data may take up to 3 minutes to display in Live View.

Please refer to Form 9215 v1 for instructions on how to access Live View.

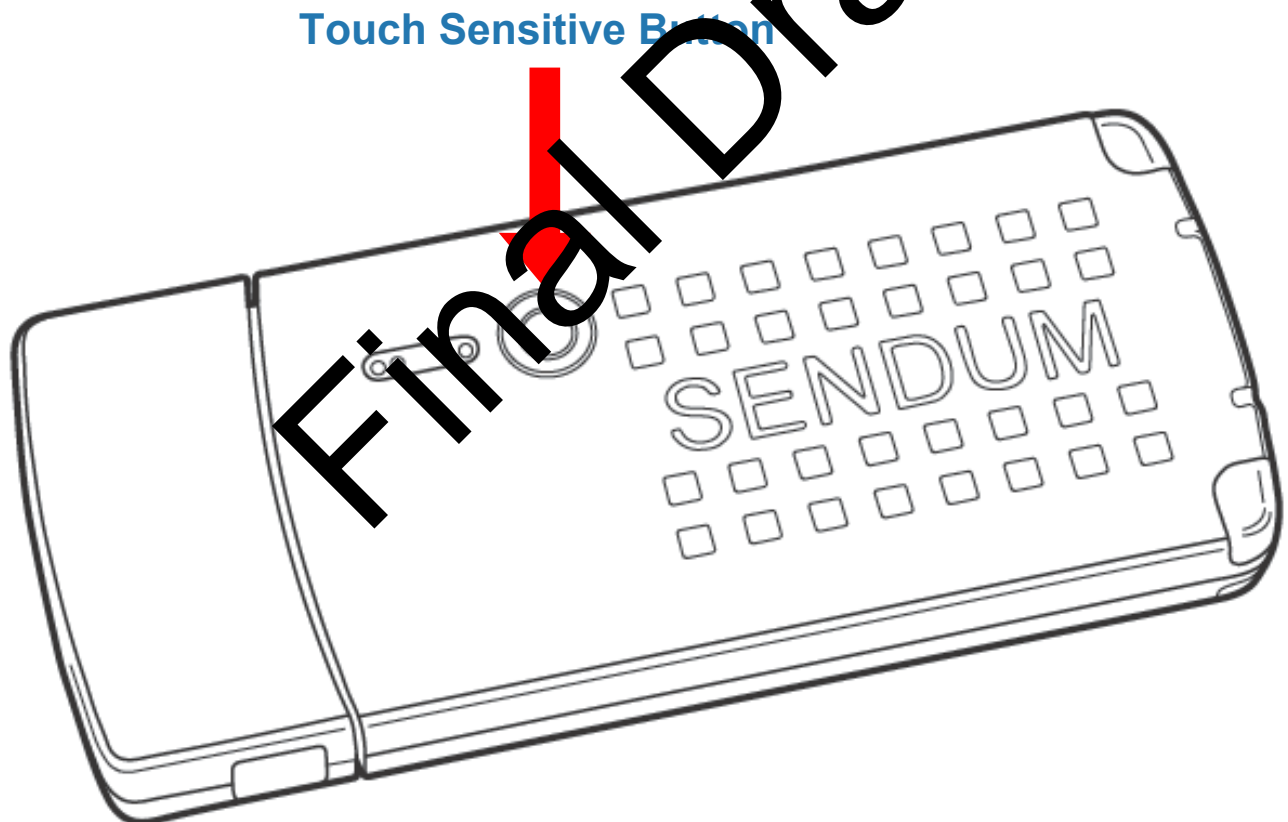

24-HOUR HELP LINE
(949) 470-2305

Cryoport, Inc. • 17305 Daimler St. • Irvine, CA 92614
Phone: 1-949-470-2300 Fax: 1-949-470-2306
www.cryoport.com • E-mail Contact: cs@cryoport.com